

PartialCAD 3.3 Chemnitz

Doing the splits





Powerful
new automatic
split denture
feature

PartialCAD 3.3 Chemnitz at a glance

PartialCAD 3.3 Chemnitz delivers fantastic new tools to simplify the design process, including a highly automated, guided split denture workflow*!


Get next-level designs with our powerful automatic split denture feature

- **Design** an implant-based suprastructure in DentalCAD® and generate a precisely fitting underlying bar **with just a few clicks** in PartialCAD
- **Easily and accurately split your All-on-X design** while also allowing manual adjustments
- **Enjoy less complexity** compared to older convoluted workflows involving tools from different vendors
- **Speed up production** and deliver high-quality results with confidence
- **Geared toward PartialCAD users or those with no extensive experience with the exocad Bar Module**, ideal for labs desiring high productivity through standardized designs and automated workflows

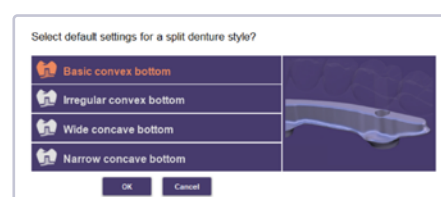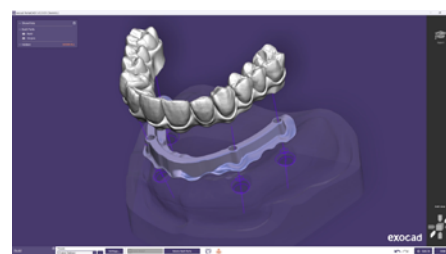

Flatter bar tops for a superior denture fit

- **Configure** the top surface for a **flawlessly level result**
- Shape and height of the bar plane are controlled by curve points attached to it, which **define the contour and elevation of the surface**

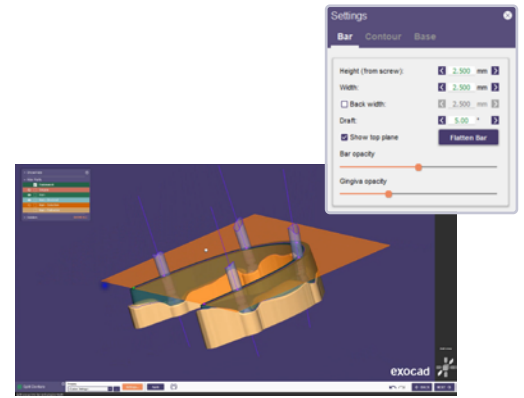


Accelerate your workflow with the enhanced post tool

- **Subtract material** across multiple parts **simultaneously** and add posts
- Get more control over design refinements and **reduce the need for repetitive changes**




Advanced visualization tools for more control

- Detect areas in the gingiva design that are too thin and automatically add material to **ensure mechanical strength**
- Save time with automation instead of manual adjustments and **feel confident about the minimal thickness settings**

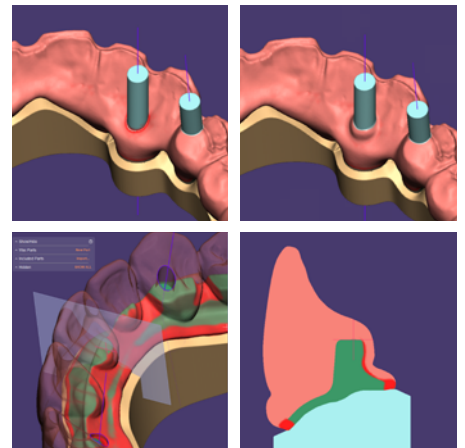

Manage your preferences like a pro

- **Save your split denture settings** so you can reuse them across cases
- Make your work **more consistent and efficient**

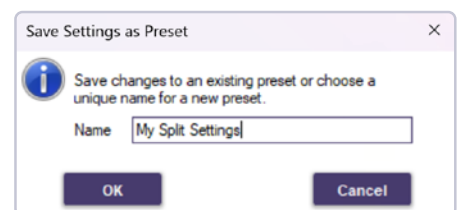

[exocad.com/partialcad-chemnitz](https://www.exocad.com/partialcad-chemnitz)

Find your exocad reseller here: [exocad.com/resellers](https://www.exocad.com/resellers)

exocad GmbH

Some products may not be regulatory cleared/released for sale in all markets. Please contact your local exocad reseller for current product assortment and availability. Territorial coverage of the registered trademarks is available on our webpage: [exocad.com/imprint](https://www.exocad.com/imprint)